


PRODUCT INFORMATION

VEGF121 (Vascular endothelial growth factor 121), Mouse

Catalog number	C02071-5UG / C2071-20UG / C02071-100UG
Package	5 µg / 20 µg / 100 µg
Description	Vascular endothelial growth factor (VEGF), originally known as vascular permeability factor (VPF), is a signal protein produced by cells that stimulates the formation of blood vessels. VEGF is required during embryogenesis to regulate the proliferation, migration, and survival of endothelial cells. In adults, VEGF functions mainly in wound healing and the female reproductive cycle. Pathologically, it is involved in tumor angiogenesis and vascular leakage. Circulating VEGF levels correlate with disease activity in autoimmune diseases such as rheumatoid arthritis, multiple sclerosis and systemic lupus erythematosus. VEGF is induced by hypoxia and cytokines such as IL-1, IL-6, IL-8, oncostatin M and TNF-alpha.
Source	<i>Escherichia coli</i>
Sequence	MAPTTEGEQKSHEVIKFM DVYQRSYCRPIETLVDIFQEYPDEIEYIFKPSCVPLM RCAGCCNDEALECVPTSES NITMQIMRIKPHQSQHIGEMSFLQHSRCECRPKKD RTKPEKCDKPRR with polyhistidine tag at the C-terminus
Endotoxin level	<0.1 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to induce proliferation in HUVEC cells. The ED ₅₀ for this effect is <3 ng/mL.
Purity	>95% as determined by SDS-PAGE. Ni-NTA chromatography
Formulation	The protein was lyophilized from a solution containing 1X PBS, pH 8.0.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Storage	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
Note	Please use within one month after protein reconstitution.

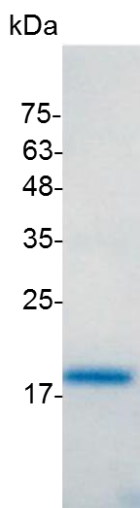

SDS-PAGE analysis of recombinant mouse VEGF121

For Research Use Only.