


Human ACE2 (Angiotensin-converting enzyme 2), Fc tag, HEK293

v. 230201

Catalog number	C11015-5UG / C11015-20UG / C11015-100UG
Package	5 µg / 20 µg / 100 µg
Source	HEK293
Sequence	<p> QSTIEEQAKTFLDKFNHEAEDLFYQSSLASWNYNTNITEENVQNMNNAGDKWS AFLKEQSTLAQMYPLQEIQNLTVKQLQALQQNGSSVLSKSKRLNTILNTMS TIYSTGKVCNPDNPQECLLLEPGLNEIMANSLDYNERLWAWESWRSEVKGKQLR PLYEEYVVLKNEMARANHYEDYGDYWRGDYEVNGVDGYDYSRGLIEDVEHT FEEIKPLYEHLHAYVRAKLMNAYPSYISPIGCLPAHLLGDMWGRFWTNLYSLTV PFGQKPNIDVTDAMVDQAWDAQRIFKEAEKFFVSVGLPNMTQGFWENSMLTD PGNVQKAVCHPTAWDLGKGDFRILMCTKVTMDDFLTAHHEMIGHIQYDMAYAA QPFLLRNGANEGFHEAVGEIMSLSAATPKHLKSIGLLSPDFQEDNETEINFLKQ ALTIVGTLPFTYMLEKWRWMVFKGEIPKDQWMKKWWEMKREIVGVVEPVPHD ETYCDPASLFHVSNDYSFIRYYTRTLYQFQFQEALCQAAKHEGPLHKCDISNST EAGQKLFNMLRLGKSEPWTLALENVVGAKNMNVRPLLNYFEPLFTWLKDQNK NSFVGWSTDWSPYADQSIKVRISLKSALGDKAYEWNENEMYLFRSSVAYAMR QYFLKVKNQMILFGEEDVRVANLKPRIFFVTAPKNVSDIIPRTEVEKAIRMS RSRINDAFRLNDNSLEFLGIQPTLGPPNQPPVS with human IgG1 Fc tag at the C-terminus </p>
Endotoxin level	<0.1 EU per 1 µg of the protein by the LAL method.
Purity	>98% as determined by SDS-PAGE.
Formulation	The protein was lyophilized from a solution containing 1X PBS containing, pH 7.4.
Reconstitution	It is recommended to reconstitute the lyophilized protein in sterile H ₂ O to a concentration not less than 100 µg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.
Storage	Lyophilized protein should be stored at -20°C. Upon reconstitution, protein aliquots should be stored at -20°C or -80°C.
Note	Please use within one month after protein reconstitution.

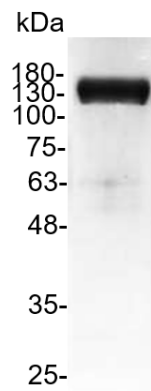

SDS-PAGE analysis of recombinant human ACE2

For Research Use Only.