

Catalog number	C01026-GMP-100 / C01026-GMP-1000
Package	100 µg / 1 mg
Description	<p>Interleukin-21 (IL-21) belongs to the IL-15/IL-21 family. It is a cytokine with immunoregulatory activity. Cytokines are proteinaceous signaling compounds that are major mediators of the immune response. They control many different cellular functions including proliferation, differentiation and cell survival/apoptosis but are also involved in several pathophysiological processes including viral infections and autoimmune diseases. IL-21 is a cytokine that has potent regulatory effects on cells of the immune system, including natural killer (NK) cells and cytotoxic T cells that can destroy virally infected or cancerous cells. This cytokine induces cell division/proliferation in its target cells.</p>
Expression System	<i>Escherichia coli</i>
Species of Origin	Human
Affinity Tag	His Tag (C-term)
Sequence	Gln32-Ser162
Endotoxin level	<0.05 EU per 1 µg of the protein by the LAL method.
Activity	Measure by its ability to enhance IFN gamma secretion in NK-92 cells. The ED ₅₀ for this effect is <10 ng/mL.
Purity	>95% as determined by SDS-PAGE analysis.
Mycoplasma	Not detected
Form	Lyophilized
Storage Buffer	Lyophilized from a 0.2 µm filtered solution of PBS, pH 8.0.
Reconstitution	<p>It is recommended to reconstitute the lyophilized protein in sterile H₂O to a concentration not less than 0.5 mg/mL and incubate the stock solution for at least 20 min to ensure sufficient re-dissolved.</p>
Stability & Storage	<p>This product is stable after storage at:</p> <ul style="list-style-type: none">• -20°C for 12 months in lyophilized state from date of receipt.• -20°C or -80°C for 1 month under sterile conditions after reconstitution. <p>Avoid repeated freeze/thaw cycles.</p>


Croyez GMP® recombinant proteins are manufactured in ISO 13485:2016 and GMP-certified facility. The processes include:

Specification

- Animal-free reagent and laboratory
- Manufactured and tested under GMP guideline
- Testing and traceability of raw material
- Records of the maintenance and equipment calibration
- Personnel training records
- Batch-to-batch consistency
- Documentation of QA control and process changes
- Manufactured and tested under an ISO 13485:2016 certified quality management system
- Stability monitor of product shelf-life

Reference

1. Spolski R and Leonard W.J. (2014) Nature Reviews Drug Discovery 13(5):379-95.
2. Leonard WJ, Spolski R. (2005) Nat Rev Immunol. 5,9: 688-98.
3. Brandt K. et al. (2007) Cytokine Growth Factor Rev. 18,3-4: 223-32.
4. Sim GC, Radvanyi L. (2014) Cytokine Growth Factor Rev. 25,4: 377-90.

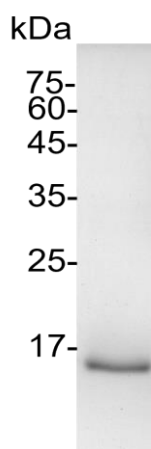

SDS-PAGE analysis of Croyez GMP® IL-21, Human

For Research Use Only. Not for use in diagnostic or therapeutic procedures.